

PRODUCT SPECIFICATIONS

Syringe Size	1 ml - 20 ml single-use syringe
Infusion Rate	1 mm/hr - 99 mm/hr
Alarms	Battery Depleted, Mechanical Error, Occlusion Syringe Disengaged, Completion.
Alerts	Battery Low, User Prompt, Syringe Not Installed, Near End
Accuracy	± 2 %
Bolus Volume	Programmable from 0.1 ml to the capacity of the syringe size.
PCA Function	The PCA volume is programmable between 0.1 ml to 10 ml in 0.1 ml increments; the locking interval is programmable between 5 minutes to 60 minutes in 5-minute increments.
OLED Display	On-screen messages of: syringe size; rate is set by mm/hr and displayed in ml/hr simultaneously; alarm messages, running status.
"LOW BATT" LED Indicator	Flashing yellow/red indicator when battery is low or depleted
Maximum Infusion Pressure	≤ 0.12 MPa (with 10 ml syringe)
Occlusion Alarm Pressure	0.12 MPa (with 10 ml syringe)
Safety Approvals	CE 0086, IEC 60601 - 1 - 8: 2003, Type CF, IPX1
Battery Size	AA x 3
Dimensions	166(W) x 60(D) x 30(H) mm
Weight	≤ 180g (battery weight not included)

canafusion 

Canafusion Technologies Inc.

301-2999 Underhill Ave. Burnaby, BC Canada. V5A 3C2
Tel: +1(604) 568.0936 / +1(604) 568.6977 Fax: +1(604) 677.5295

info@canafusion.com

www.canafusion.com


canafusion 



CA-500 Ambulatory Syringe Pump


CE
0086

CA-500 Ambulatory Syringe Pump

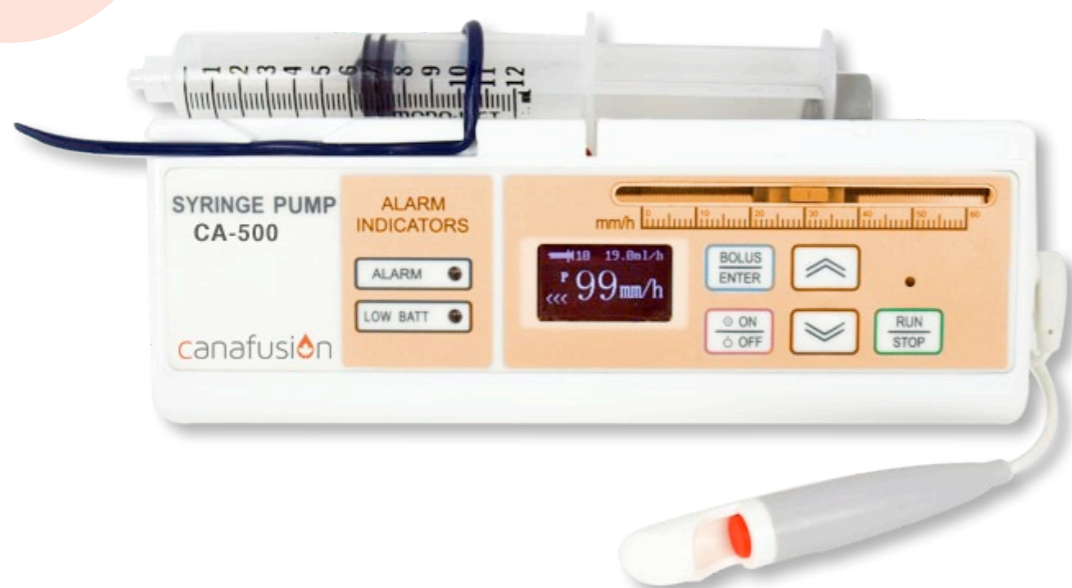

Our CA-500 ambulatory syringe pump is an accurate, reliable infusion solution for Thalassemia (desferral) treatment, pain management, neonatal care and many other ambulatory treatments with low rate small volume. It's compatible with any brand syringe from 1 to 20ml and is programmable in mm/h. The pump can also automatically convert & display the ml/hr rate based on mm/h rate of the calibrated syringe.

SAFETY & ACCESSORIES

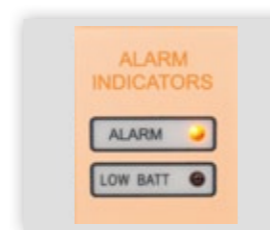


Alarm Indicators
An immediate identification of pump's alarm and alert messages.



Bolus Cable (PCA)
Patients can click bolus button to inject extra dose for quick relief of pain.



Bright OLED Screen
Bright & Energy saving OLED screen displays both infusion rate and alarm message.

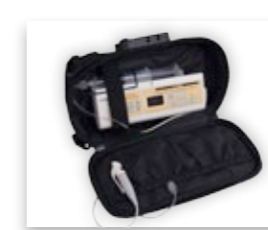

Carrying Bag
The carrying bag comes with programmable lock and exit for the PCA cable and IV line for patient.



Syringe Cover
Protects the syringe during the infusion to guarantee a safe & accurate infusion.



Package Contents
The package includes: a carrying bag with programmable lock for ambulatory use, syringe covers, operating instructions.

Great for Home Care or Ambulatory Treatment

APPLICATIONS



Thalassemia



Parkinson (Apomorphine hydrochloride)



Neonatal Care



Immunity Deficiency (Globulin)

EASE OF USE

- Lightweight and compact.
- Digital programming with a bright OLED screen.
- Clear display of current infusion rate and alarm message.
- Powered by 3 X AA batteries.

